

Amount of tooth movement achieved at the end of the treatment

 **Blue: Post-Treatment**




 **White: Pre-Treatment**

V Hello VISHWAKRIT CHORADIA

Case ID: 2401071894

Get ready for the memorable journey towards your perfect smile!! Your customised treatment plan is ready based on your desire. Let's explore the digital design of your comprehensive treatment plan.

PRE-TREATMENT VIEW

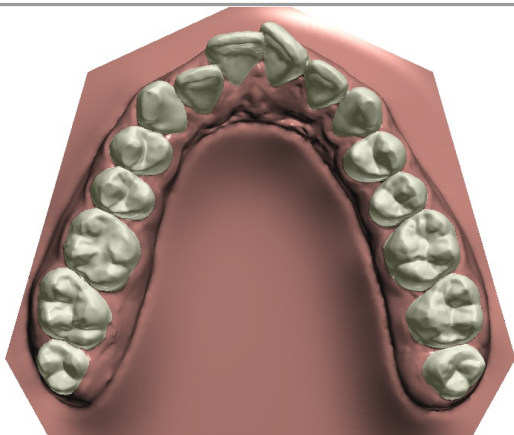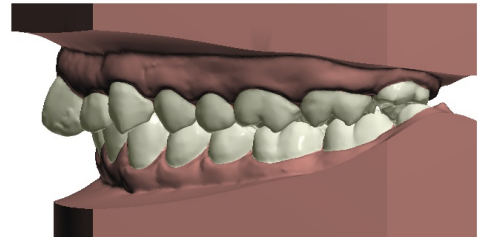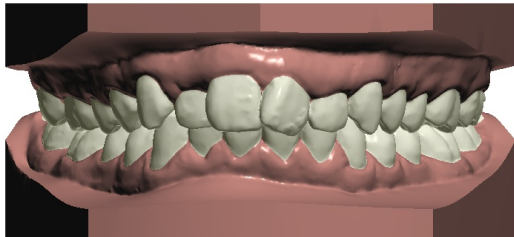

ANNEXURE-1

Tooth width analysis

11	7.66 mm	21	7.56 mm	Actual arc length: 113.44 mm
12	6.22 mm	22	5.75 mm	
13	6.34 mm	23	6.57 mm	
14	6.07 mm	24	5.90 mm	
15	6.25 mm	25	5.27 mm	
16	9.30 mm	26	9.14 mm	
17	8.95 mm	27	9.06 mm	
18	6.01 mm	28	7.38 mm	
48	5.01 mm	38	4.53 mm	Actual arc length: 102.39 mm
47	8.18 mm	37	8.98 mm	
46	8.96 mm	36	9.57 mm	
45	5.78 mm	35	6.69 mm	
44	6.10 mm	34	6.82 mm	
43	5.95 mm	33	5.73 mm	
42	4.98 mm	32	5.39 mm	
41	4.64 mm	31	5.09 mm	

Space Analysis

MAXILLARY SPACE ANALYSIS

- Total space available: 63.60
- Sum of incisor widths: 27.20

Moyers

- Sum of left canine & premolars: 20.96
- Sum of right canine & premolars: 20.96
- Total space required: 69.11
- Discrepancy: -5.52

Tanaka & Johnston

- Sum of left canine & premolars: 21.05
- Sum of right canine & premolars: 21.05
- Total space required: 69.30
- Discrepancy: -5.70

MANDIBULAR SPACE ANALYSIS

- Total space available: 57.17
- Sum of incisor widths: 20.09

Moyers

- Sum of left canine & premolars: 20.46
- Sum of right canine & premolars: 20.46
- Total space required: 61.01
- Discrepancy: -3.84

Tanaka & Johnston

- Sum of left canine & premolars: 20.55
- Sum of right canine & premolars: 20.55
- Total space required: 61.19
- Discrepancy: -4.02

Bolton analysis

ANTERIOR BOLTON ANALYSIS

- Maxillary 6 width: 40.11 mm
- Mandibular 6 width: 31.78 mm
- Ratio: 0.79

Maxill. 6	Mandib. 6	Ratio
40.00	30.90	0.77
40.50	31.30	0.77
41.00	31.70	0.77
41.50	32.00	0.77
42.00	32.40	0.77
42.50	32.80	0.77
43.00	33.20	0.77
43.50	33.60	0.77
44.00	34.00	0.77
44.50	34.40	0.77
45.00	34.70	0.77
45.50	35.10	0.77
46.00	35.50	0.77
46.50	35.90	0.77
47.00	36.30	0.77
47.50	36.70	0.77
48.00	37.10	0.77
48.50	37.40	0.77
49.00	37.80	0.77
49.50	38.20	0.77
50.00	38.60	0.77
50.50	39.00	0.77
51.00	39.40	0.77
51.50	39.80	0.77
52.00	40.10	0.77
52.50	40.50	0.77
53.00	40.90	0.77
53.50	41.30	0.77
54.00	41.70	0.77
54.50	42.10	0.77
55.00	42.50	0.77


OVERALL BOLTON ANALYSIS

- Maxillary 12 width: 82.04 mm
- Mandibular 12 width: 75.70 mm
- Ratio: 0.92


Maxill. 12	Mandib. 12	Ratio
85.00	77.60	0.91
86.00	78.50	0.91
87.00	79.40	0.91
88.00	80.30	0.91
89.00	81.30	0.91
90.00	82.10	0.91
91.00	83.10	0.91
92.00	84.00	0.91
93.00	84.90	0.91
94.00	85.80	0.91
95.00	86.70	0.91
96.00	87.60	0.91
97.00	88.60	0.91
98.00	89.50	0.91
99.00	90.40	0.91
100.00	91.30	0.91
101.00	92.20	0.91
102.00	93.10	0.91
103.00	94.00	0.91
104.00	95.00	0.91
105.00	95.90	0.91
106.00	96.80	0.91
107.00	97.80	0.91
108.00	98.60	0.91
109.00	99.50	0.91
110.00	100.40	0.91

Front View

Pre - Treatment

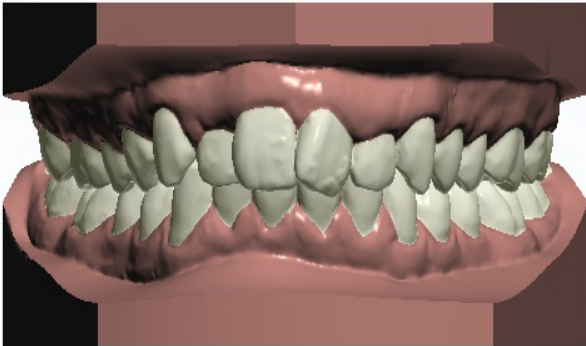


Post - Treatment

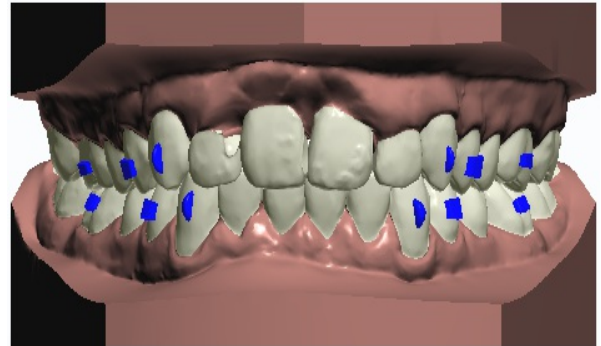


Right View

Pre - Treatment

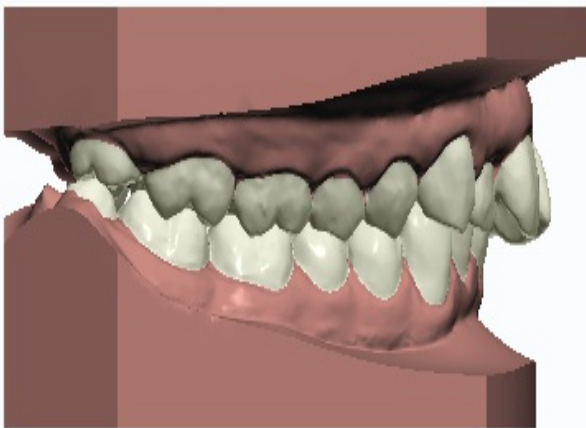


Post - Treatment

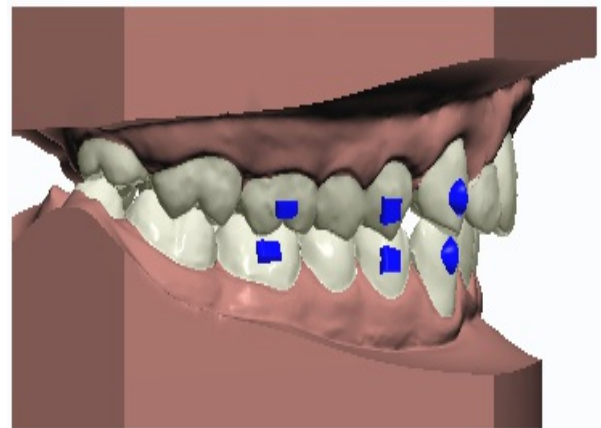


Left View

Pre - Treatment

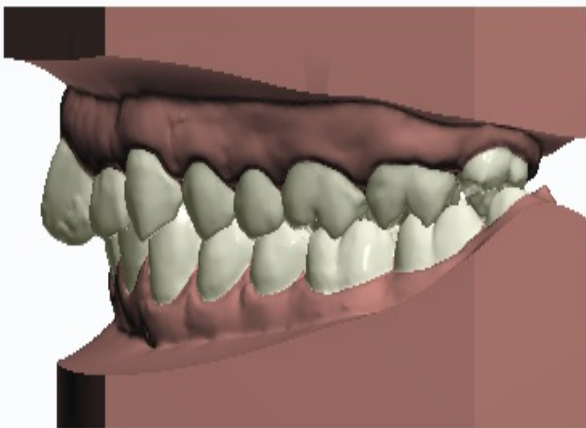


Post - Treatment

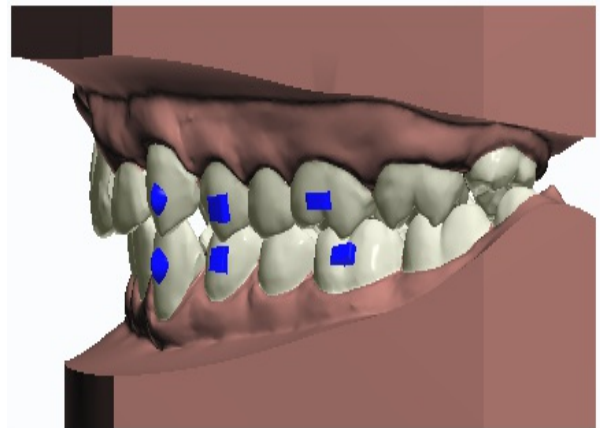


Upper View

Pre - Treatment

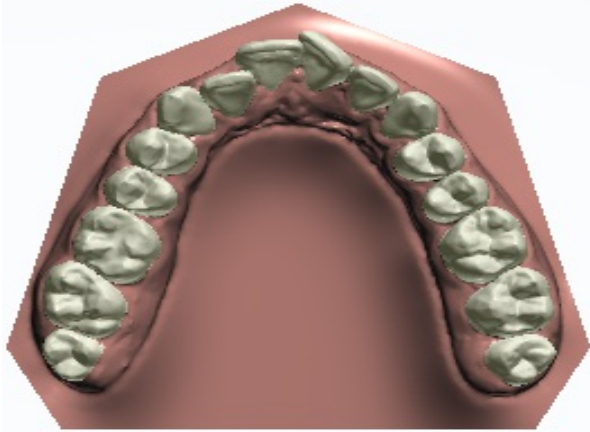

Post - Treatment

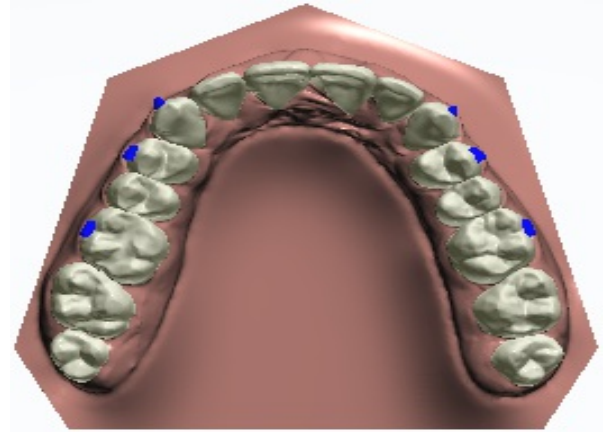


Lower View

Pre - Treatment



Post - Treatment

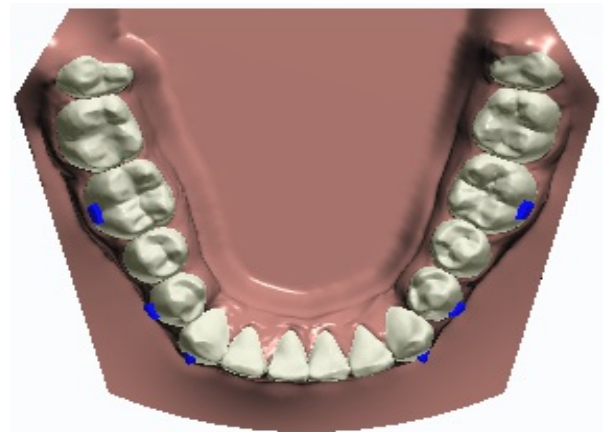

ANNEXURE-3

Teeth Movement Report

upload plan

Tooth number	Rotation	Angulation	Inclination	Left/ Right	Extrusion/ Intrusion	Forward/ Backward	IPR Distal	IPR Mesial	Long Axis	Circumference
16							Manual	Manual	82.2 deg	37.1 mm
15			2.0 deg	-0.9 mm			Manual	Manual	81.6 deg	26.9 mm
14			4.0 deg	-1.7 mm			Manual	Manual	77.5 deg	26.2 mm
13	-15.0 deg			-2.1 mm		0.9 mm	Manual	Manual	76.8 deg	23.8 mm
12	9.2 deg		1.0 deg	-2.4 mm		-0.1 mm			75.3 deg	20.1 mm
11	-1.0 deg		-5.0 deg	-0.4 mm	1.4 mm	-2.0 mm			76.0 deg	23.6 mm
21	-38.3 deg			-3.6 mm	1.8 mm	-1.1 mm	Manual	Manual	68.1 deg	23.3 mm
22	-10.0 deg			-2.6 mm	0.5 mm	0.4 mm	Manual	Manual	74.0 deg	19.2 mm
23	-8.0 deg			-2.2 mm	0.3 mm	0.8 mm	Manual	Manual	82.3 deg	23.0 mm
24				-1.5 mm		0.8 mm	Manual	Manual	86.7 deg	26.4 mm
25				-0.8 mm		0.4 mm	Manual	Manual	84.2 deg	25.8 mm
26							Manual	Manual	77.5 deg	36.4 mm
33	-16.0 deg		3.0 deg	0.5 mm	1.1 mm	0.3 mm			80.4 deg	21.2 mm
32	-10.0 deg			-0.1 mm	0.9 mm	-0.5 mm			72.0 deg	17.3 mm
31	-4.0 deg		-4.0 deg	0.3 mm	1.4 mm	-0.4 mm	Manual	Manual	68.0 deg	16.2 mm
41	-1.0 deg		-5.0 deg	0.2 mm	1.5 mm	-0.2 mm	Manual	Manual	68.4 deg	15.9 mm
42	-11.0 deg			-0.5 mm	0.9 mm	0.1 mm	Manual	Manual	70.3 deg	16.1 mm
43	-13.0 deg			0.1 mm	1.0 mm	0.3 mm	Manual	Manual	82.3 deg	22.0 mm
44						0.1 mm			79.0 deg	23.5 mm
45			1.0 deg						83.6 deg	24.9 mm

Do reach out to our team if you have any further questions and we'd be happy to have our executives speak to you. Have a great day, and remember, keep smiling. 😊

Best Regards
(Team Rejove32 Aligners)